



MEDICAL-X



# Trauma leg

Task Trainer

# Trauma Leg Task Trainer

Practice first aid measures and manage a major hemorrhage with the Trauma Leg.

The trainer represents an injured leg with laceration, bullet entry, and bullet exit wounds.

The bleeding system enables the flow of blood from the wounds.

You can apply direct pressure, pack the wound with gauze, place a pressure dressing, and apply the tourniquet to really stop the bleeding.

Made with a soft silicone coating, the Trauma Leg reproduces the feeling of touching the human skin.





## What you need



**Simulated blood**

*INCLUDED*



**Application bottle**

*INCLUDED*



**Tourniquet**

*NOT INCLUDED*

## Preparation



**Step 1:** Insert the bottle nozzle into one of the tubes at the top side of the Trauma Leg. Each tube is connected with one of the three wounds.

**Step 2:** Fill the bottle with simulated blood.

**Step 3:** Attach the application bottle to the

Trauma Leg Trainer by gently inserting the nozzle of the bottle into one of the tubes at the top side of the wound packing trainer. Each tube is connected with one of the three wounds.

**Step 4:** Apply some force to make sure the nozzle is secured tightly.

**Step 5:** Squeeze the application bottle to push the simulated blood through the tubes into the wound.



Detachable tubes



## Bleeding

The bleeding system can be activated by inserting simulated blood into the tube.


**Each tube is connected to a specific wound.** See the next image for the wounds with their corresponding tubes.

### Tip:

Use our new automatic pump

**Hemorrhage-X** for an easier and more realistic bleeding experience.

Contact us to know more about!

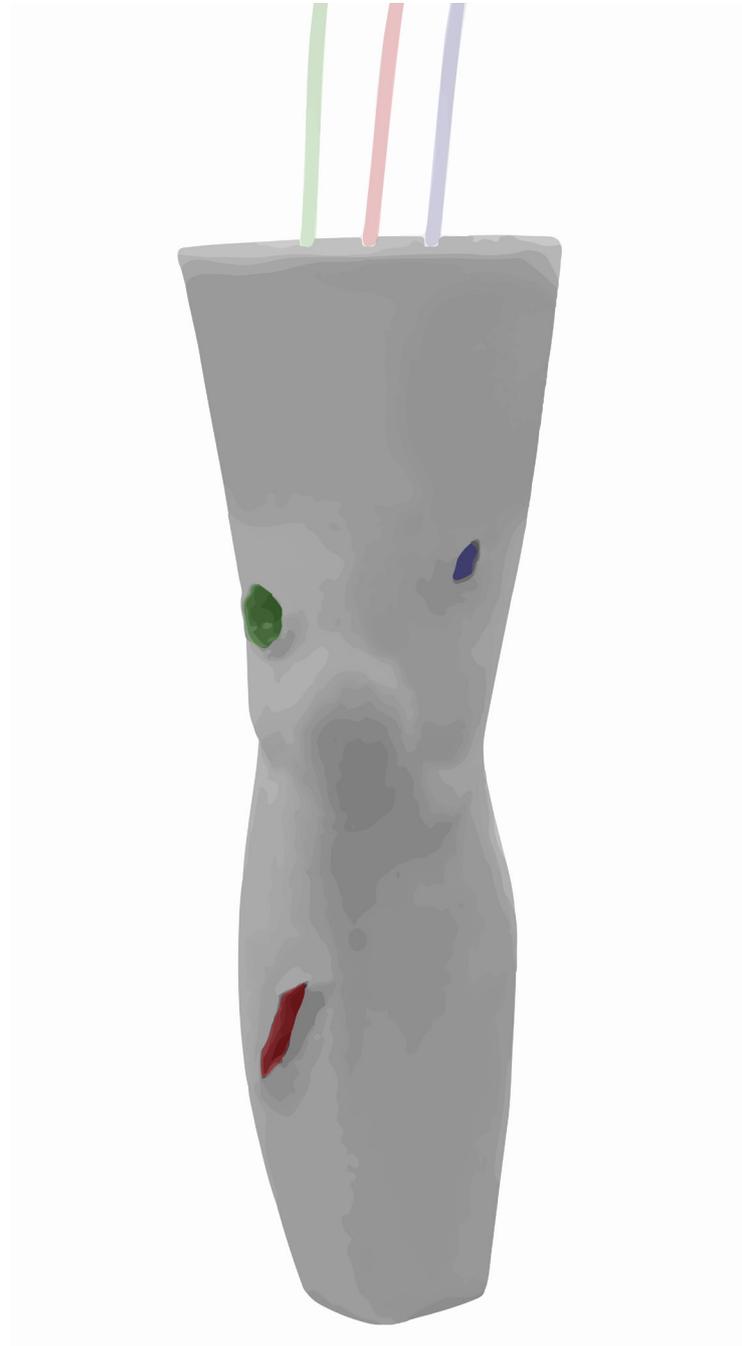

# Cleaning



**Step 1:** Detach the application bottle from the tube and pour back the remaining content into the original artificial blood bottle.

**Step 2:** Fill the application bottle with a light soap solution or mild domestic cleaners without alcohol.

**Step 3:** Provide a container or sink to dispose of the remaining blood from the wounds. Empty the wounds by holding the leg sideways and upside-down.

**Step 4:** Reattach the application bottle and flush the tubes and wounds with an abundant amount of water.

**Step 5:** Use a soft cloth to clean the wounds and the skin.

**Step 6:** Repeat these steps until the Trauma Leg Task trainer is properly cleaned.

**Step 7:** Make sure the inside of the tubes and wounds are as good as dry. A little wetness will dry by itself when the Trauma Leg Task Trainer is left in open air.

**Step 8:** Dry off the outside of the Trauma Leg Task Trainer.

**Step 9:** Detach the tubes by pressing in the grey ring, this will release the pin and it can be pulled out.

# Storage

Store the Trauma Leg Task Trainer in a dry environment at room temperature, out of direct sunlight.



# Features List

## Dimensions

- 48 cm (18.9 inches) in length
- Diameter 18 cm (7 inches)

## Specifications

- The bleeding system enables the flow of blood from the wounds.
- Tourniquet application stops the bleeding.

## Materials

Silicone and plastic representing soft tissues.



## CONTACT US

Contact us for more information about our product and the ordering process.

🌐 [www.medical-x.com](http://www.medical-x.com)

✉ [info@medical-x.com](mailto:info@medical-x.com)

☎ +31 10 341 5927

📍 Schiedamsedijk 77, 3011 EM  
Rotterdam, The Netherlands